

Premium testing solutions for your health journey

convenient at-home sampling and
precise clinical laboratory results

PRODUCT CATALOG

2026

Cortisol / Omega Ratio / Vitamin D / NAD Level / HbA1c / Testosterone



Get Your
Results

Ship to
Our Lab

Collect Your
Sample

Order
Your Kit





TABLE OF CONTENTS

About Us	4
Development of Laboratory Methods	7
CORTISOL TEST	8
OMEGA RATIO TEST	12
VITAMIN D TEST	16
NAD LEVEL TEST	20
HbA1c TEST	24
TESTOSTERONE+ TEST	28

About Us

Your health is our priority

At Bloodo, we're more than just a health company. We're passionate about giving you the power to understand your health on a deeper level. Our convenient at-home blood sample collection kits, like our NAD Level Test, Omega Ratio Test, Vitamin D and HbA1c tests, put advanced technology right in your hands. We don't just send results; we provide insights at the cellular level, empowering you to make informed choices about your well-being.

We believe that understanding your health shouldn't be complicated. That's why we're with you every step of the way, offering resources and knowledge to guide you on your journey. Your health is your most valuable asset – let us help you protect it.

Advanced health insights from the comfort of home

Bloodo revolutionizes personal health monitoring with our innovative at-home blood testing solutions. We utilize Dried Blood Spot (DBS) technology, a cutting-edge approach that offers several advantages:

- **Molecular-level insights:**

DBS enables the detection of key health markers at the molecular level, providing in-depth understanding of your well-being.

- **Convenience and simplicity:**

DBS requires only a small blood sample, easily collected and shipped to our laboratory for analysis using advanced techniques like LC-MS/MS.

- **Enhanced stability:**

Dried blood samples eliminate the need for specialized storage, simplifying the process for both you and our laboratory.

- **Safety:**

DBS minimizes the risk of contamination during infectious disease diagnosis, ensuring safety for all involved.





Innovation and precision: welcome to LifeLab1

At Bloodo we're committed to making health monitoring more accessible, accurate, and safe, empowering you to make informed decisions about your well-being.

Nestled in Vilnius, Lithuania, LifeLab1 is a state-of-the-art clinical laboratory where cutting-edge technology meets scientific expertise. Utilizing advanced instruments like liquid chromatography/mass spectrometry (LC/MS), gas chromatography/mass spectrometry (GC/MS), and high-performance liquid chromatography (HPLC) systems, along with automated sample preparation, LifeLab1 delivers highly accurate and reliable health insights.

Fully certified and licensed for the provided services: NM-334, NR-4864, adhering to ISO 15189 requirements.

Methods: colorimetric assay, ELISA (enzyme-linked immunosorbent assay), LC-MS/MS (liquid chromatography – tandem mass spectrometry); GC/MS (gas chromatography – mass spectrometry) ; HPLC (high performance liquid chromatography).

Powered by experts, driven by quality


LifeLab1 team of dedicated scientists, researchers, medical biologists and managers ensures that your results are both comprehensive and actionable. Laboratory adheres to ISO 15189 standards, providing you with the utmost confidence in the accuracy of your test results.

LifeLab1 is committed to excellence through:

- External quality control programs
- A professional and highly qualified team
- State-of-the-art equipment
- Rigorous quality procedures
- Collaboration with leading scientists



DEVELOPMENT OF LABORATORY METHODS

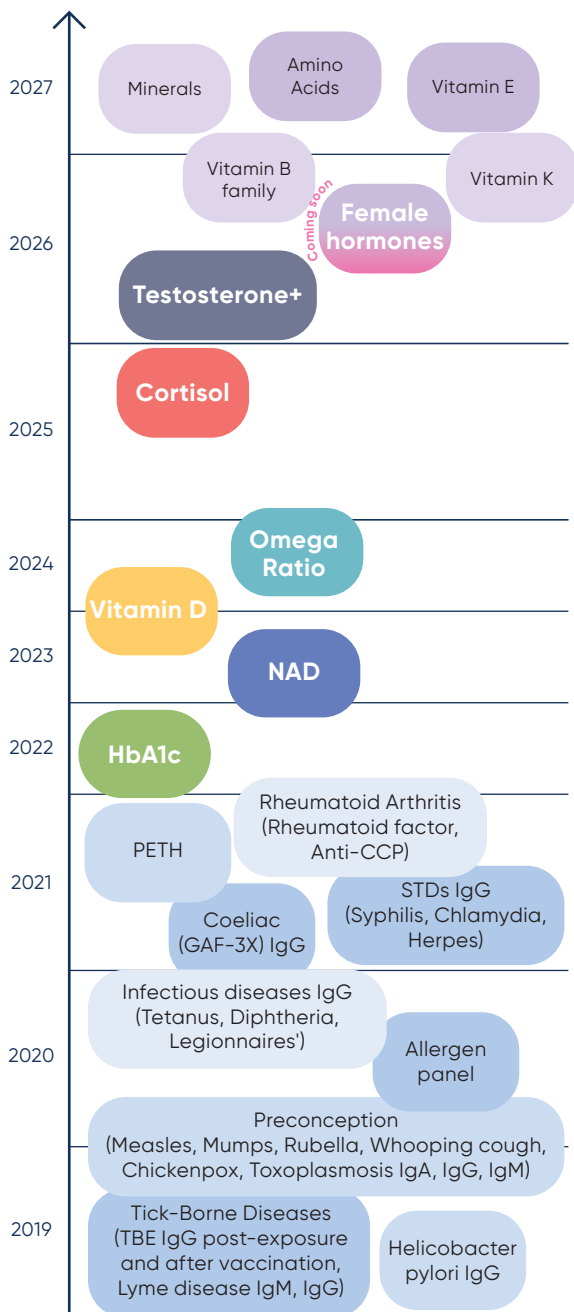

CORTISOL TEST

Uncover and balance your stress pattern



Stress
response



Energy
regulation



Sleep-wake
cycle



Metabolic
balance

Order your
test here:



- Monitor stress adaptation and resilience
- Support energy and recovery balance
- Evaluate sleep-wake cycle regulation
- Detailed analysis: morning and evening cortisol levels and their balance.

What is Cortisol?

Cortisol is a glucocorticoid hormone produced in the adrenal cortex from cholesterol. Although often referred to as the “stress hormone”, it plays a key role in maintaining physiological balance by regulating metabolism, immune responses, sleep–wake cycles, and brain function. Cortisol secretion is controlled by the hypothalamic–pituitary–adrenal (HPA) axis and follows a circadian rhythm, peaking after waking and declining throughout the day. Disruption of this rhythm, rather than cortisol itself, is associated with adverse health effects.

Why test Cortisol?

Cortisol testing provides information about adrenal function, the activity of the hypothalamic–pituitary–adrenal (HPA) axis, and the body’s response to stress. Since cortisol follows a strong daily rhythm, assessing its morning peak and evening decline is particularly valuable for detecting physiological imbalances. Short-term increases in cortisol are a normal and adaptive response. However, chronically elevated or insufficient cortisol levels, or a dysregulated daily rhythm, can contribute to physiological stress and increase the risk of organ system dysfunction.

Cortisol testing applications:

- Assessment of pituitary or adrenal function activity.
- Insight into stress-related conditions, including burnout, sleep disturbances, chronic fatigue, and altered stress reactivity.
- Monitoring treatment/intervention responses, such as glucocorticoid therapy, stress-management interventions, lifestyle changes, or circadian therapies, ensuring normalized peaks/nadir.
- Research applications, where cortisol serves as a biomarker of stress resilience, circadian health, metabolic adaptability, or immune function in various populations.

How does our test work?

With just a few drops of blood, we can accurately determine your cortisol levels and assess how well your body is managing daily stress. Our test is performed using enzyme-linked immunosorbent assay (ELISA) – a highly sensitive and precise analytical method widely used in clinical diagnostics. This approach allows for reliable cortisol measurement even from dried blood spot (DBS) samples, making the entire process simple and accessible.

Our at-home method captures cortisol at its natural extremes – in the morning, shortly after waking, when levels are highest, and in the evening, just before sleep, when they are lowest. This allows us to assess the true daily cortisol rhythm and calculate a more accurate morning-to-evening ratio – an insight that standard testing workflows frequently overlook.

Specifications of *Bloodo™* Cortisol Test

Parameter	Specification
Method principle	Enzyme-Linked Immunosorbent Assay (ELISA) for quantitative cortisol measurement from dried blood spot (DBS) samples.
Suitability	For human dried blood spot analysis.
Linearity	≥ 0.99
Limit of detection, whole human blood	11.11 nmol/l
Limit of quantification, whole human blood	18.52 nmol/l
Linear range, whole human blood	11.04–966 nmol/l
Shipping conditions (DBS test kits)	Ship at ambient temperature.
Storage conditions (DBS test kits)	Store at ambient temperature.
Shipping conditions (DBS Test Cards, blood sample placed)	Ship at ambient temperature. Post immediately. At the latest, on the following day. Sample stability limit – 34 days. Maximum sample shipping duration – 27 days.
Home sampling conditions	Collect two blood spot samples twice per day – one in the morning and one in the evening – to capture your natural cortisol rhythm. Morning sample: collect within 1 hour after waking. Evening sample: collect just before going to sleep. Avoid eating, caffeine, alcohol, and nicotine for at least 30–60 minutes before each sample. Allow the DBS Test Card to dry completely for a minimum of 3 hours before shipping it to the laboratory.

Test details

1. How are Cortisol results interpreted?

Healthy individuals exhibit a steep morning-to-evening slope, with cortisol levels rising post-awakening and then falling progressively. Flatter slopes indicate dysregulation, common in chronic stress, shift work, physical activity, or mood disorders.

An AM/PM ratio above 2:1 generally reflects a healthy cortisol rhythm.

A low AM/PM ratio indicates a flattened diurnal curve, where morning levels do not sufficiently exceed evening levels, signaling a disrupted daily pattern. Such profiles are associated with fatigue, cognitive difficulties, mood disturbances, poorer sleep quality, and elevated cardio-metabolic risk. A persistently flattened rhythm also correlates with reduced physical functioning and a higher burden of psychiatric symptoms.

Technical details

1. Calibration

Cortisol is quantified using a calibration curve prepared from cortisol standards added to real human whole blood. A typical calibration is performed across the concentration range of 11.04–966 nmol/L, enabling accurate determination of cortisol levels in DBS samples within this physiological and clinical span.

2. Accuracy

To evaluate analytical accuracy, cortisol levels measured in DBS samples (x-axis) were compared with corresponding serum samples collected at the same time. A subset of donors provided venous serum samples for this comparison.

The correlation between DBS and serum cortisol measurements was strong, yielding $R^2 = 0.9823$, demonstrating excellent agreement between the two matrices (Figure 1).

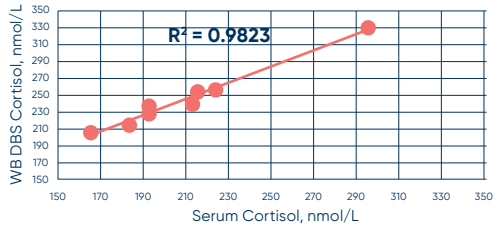


Figure 1. Comparison of cortisol concentrations measured from dried blood spot samples versus matched serum samples

3. Repeatability

Between-run repeatability of the cortisol assay was assessed by analyzing each of ten samples ten times within the same day. For samples with low cortisol concentrations, the coefficient of variation (CV) was 10.16%, whereas for samples with high cortisol concentrations, the CV was 10.59%. These results confirm that the method delivers consistent and reliable measurements across physiologically relevant cortisol ranges.

4. Reproducibility

The reproducibility of the cortisol assay was evaluated by performing multiple independent measurements on a series of samples across different days and operators. The coefficient of variation (%CV) was 10.50% for low-concentration cortisol samples and 13.36% for high-concentration samples. These results indicate that the method exhibits robust analytical performance, demonstrating high precision and consistent measurement reliability across physiologically relevant cortisol levels.

5. Stability

To assess cortisol stability, a 34-day room-temperature stability study was conducted. The results showed that cortisol concentration remained stable over the tested period, with no meaningful degradation observed. Measured concentrations stayed within acceptable stability limits, with changes not exceeding $\pm 25\%$ of the initial value (Figure 2).

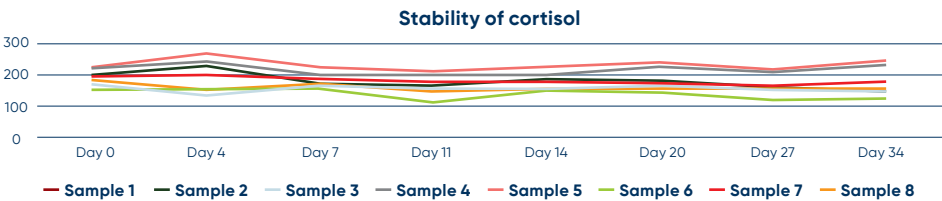


Figure 2. Stability of cortisol during 34-day storage at room temperature

OMEGA RATIO TEST

Uncover your full essential fatty acids profile



Inflammation



Cognitive
function



Heart
health

- Track cardio vitality and cognitive function
- Support inflammation management
- Check for nutritional deficiencies
- Detailed analysis: Omega 3/6/9 and others

Order your
test here:



What are Omega Fatty Acids?

Omega fatty acids are a family of polyunsaturated fats essential for many bodily functions. The three main types are omega-3, omega-6, and omega-9 fatty acids, which are important dietary fats. There are three types of omega-3s: ALA (alpha-linolenic acid) from plant sources, and EPA (eicosapentaenoic acid) and DHA (docosahexaenoic acid) from fatty fish. Omega-3s help lower triglycerides, reduce blood pressure, and decrease the risk of heart disease and stroke. Omega-6 fatty acids, found in vegetable oils, are needed in smaller amounts for growth and development. Omega-9 fatty acids, found in olive oil and nuts, are not essential but contribute to heart health and help reduce inflammation.

Why test Omega Ratio?

The Omega Ratio Test reveals potential imbalances or deficiencies in your Omega fatty acids profile. If the balance between is off, especially between omega-3 to omega-6 ratio, it can contribute to a whole bunch of health problems like heart issues, cancer risks, inflammation, autoimmune troubles, neurological conditions, and even skin issues like eczema and hair loss. Testing is a perfect way to assess your Omega fatty acid ratios and, if needed, adjust your diet or supplements for better heart health and inflammation management. If you're already taking omega-3 supplements or have made dietary changes, testing helps monitor their impact on your fatty acid balance.

What does the Omega Ratio Test measure?

Our Omega Ratio Test offers a comprehensive analysis of fatty acids within your blood, providing an in-depth view of your fatty acid balance. We analyse and report 26 individual fatty acids, along with the omega-3, omega-9, and trans fat indexes, as well as the omega-6:omega-3 and AA:EPA ratios. These key measurements help assess inflammation levels, cardiovascular risk, and overall metabolic health. This detailed approach provides a complete understanding of your personal omega profile and its potential implications for your well-being.

How does our test work?

It's amazing how much information a single drop of blood can reveal about your fatty acid profile. Accurately measuring omega fatty acid levels demands specialized laboratory expertise and techniques. Our test employs well-proven and reliable methods like gas chromatography with flame ionisation detection (GC-FID) for superior reliability, reproducibility, identification and precision. Our meticulous approach, from careful handling to superior precision, ensures results you can confidently use to guide your health choices.

Specifications of *Bloodo™ Omega Ratio Test*

Parameter	Specification	
Method principle	Gas chromatography with flame ionization detection for the analysis of fatty acids.	
Suitability	For human dried blood spot analysis.	
Repeatability, CV	Saturated	3.84 %
	Unsaturated	1.96 %
	Omega-3 index	2.93 %
	Omega-6:Omega-3	4.97 %
	AA:EPA	2.80 %
	Trans fat index	10.28 %
	Omega-9 index	5.22 %
Shipping conditions (DBS Test Kits)	For human dried blood spot analysis.	
Storage conditions (DBS Test Kits)	Store at ambient temperature.	
Shipping conditions (DBS Test Cards, blood sample placed)	Ship at ambient temperature. Post immediately. At the latest, on the following day. Sample stability limit - 42 days. Maximum sample shipping duration - 35 days.	
Home sampling conditions	Take the blood sample at a convenient time for you. Dry the DBS Test Card for 10 minutes.	

Test details

1. Purpose of the test

This test provides a comprehensive analysis of the individual fatty acid profile from a single blood sample, analyzing 26 distinct fatty acids. These include saturated fatty acids (SFAs), monounsaturated fatty acids (MUFAs), and polyunsaturated fatty acids (PUFAs). This extended fatty acid profile enables the precise calculation of several critical health biomarkers, such as the omega-3 index (normalized to red blood cell equivalent values), the omega-6:omega-3 ratio, the AA:EPA ratio, and indexes for trans fatty acids and omega-9 levels. All these parameters are calculated from a dried blood spot collected onto filter paper treated with a fatty acid-stabilizing agent.

Technical details

1. Repeatability

To assess the repeatability of our test method, samples were measured multiple times under identical conditions. The results demonstrated consistent performance across all replicates, with variations ranging from 1.15 % to 10.28 %. This range reflects the inherent variability within the measurement process, confirming that the method exhibits reliable repeatability under controlled conditions. The observed variation indicates a high level of precision, supporting the robustness of our testing approach for fatty acid composition analysis.

2. Reproducibility

The reproducibility of omega-3 index measurements was assessed by performing 6 independent measurements on 12 samples. A coefficient of variability (%CV) was calculated to be 9,97 %, reflecting the high level of precision and repeatability of our testing approach.

3. Stability

To assess the stability of omega fatty acids, a 42-day stability test was conducted. The results demonstrated that the omega-3 index remained stable at room temperature for up to 42 days, with decreases ranging from 0.9 % to 7.6 %, not exceeding the acceptable limit of 25 % (Figure 1). Other indexes, including the omega-6:omega-3 ratio, AA:EPA ratio, trans fat index, and omega-9 index, showed variations below 15 %. (Figure 1).

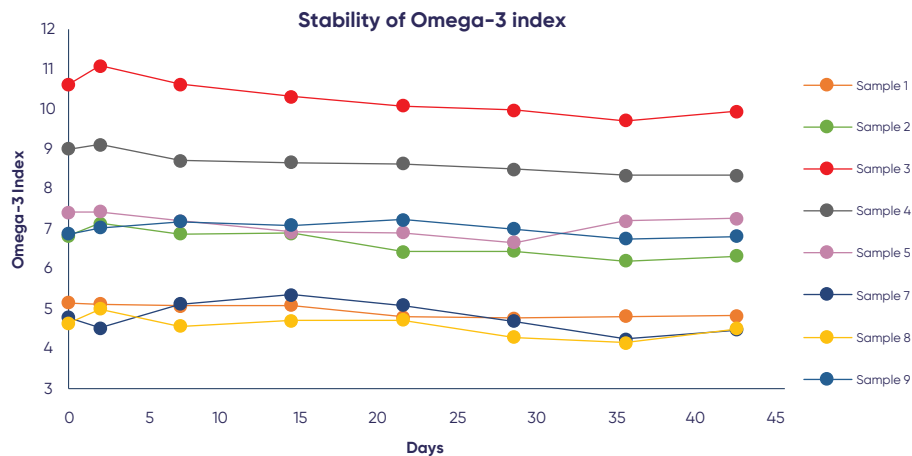


Figure 1. Stability of Omega-3 Index Over 42 Days at Room Temperature

VITAMIN D TEST

Check your sunshine levels



Immune
system



Bone and
muscle health



Mental
wellness

- Monitor immunity
- Support mental wellness
- Monitor muscle and bone strength
- Analysis of Vitamin D2 and D3

Order your
test here:



This test kit is suitable for infants. Sample collection for infants must be performed by qualified professionals.

What is Vitamin D?

Vitamin D is a fat-soluble compound essential for bone health, immune function, and overall metabolism. We call this molecule Vitamin D, but the active form, called cholecalciferol, is actually more like a hormone. It exists in two main forms: Vitamin D₃ (cholecalciferol), produced naturally by the body in response to sunlight, and Vitamin D₂ (ergocalciferol), found in some food sources. Both Vitamin D forms are with the similar structure. Vitamin D enables the body to absorb and utilise calcium to build and maintain strong bones and teeth, helping prevent rickets in children and osteomalacia in adults. Beyond bone health, it supports immune defence, regulates metabolism, and contributes to muscle strength, including the function of the heart. Insufficient Vitamin D levels are also associated with low mood and depression-like symptoms, highlighting its broad impact on physical and mental well-being.

Why test Vitamin D?

Assessing your Vitamin D levels is important, as deficiencies are widespread. Testing is especially valuable for people living in regions with low annual sunlight exposure. Not enough Vitamin D may lead to bone weakness, muscle problems, and a weakened immune system. It's important to have your Vitamin D levels tested before taking supplements. Even though it's beneficial, excessive Vitamin D can cause a dangerous buildup of calcium in your blood known as hypercalcemia. Testing helps identify if you need supplementation and guides dosage recommendations for optimal benefits.

What kind of Vitamin D do we measure?

Our Vitamin D test measures both D₃ (cholecalciferol) and D₂ (ergocalciferol) separately, ensuring precise evaluation. Our Vitamin D test measures both Vitamin D₃ (cholecalciferol) and Vitamin D₂ (ergocalciferol) separately, ensuring a precise evaluation. The report provides individual Vitamin D₃ and D₂ values, along with the total Vitamin D concentration. We measure the active form, 25-hydroxyvitamin D [25(OH)D], which is the best indicator of your body's Vitamin D stores, while excluding the less biologically active epimeric forms to ensure accuracy. We measure 25-hydroxyvitamin D [25(OH)D], the best indicator for your body's Vitamin D stores.

How does our test work?

With just a few drops of blood, we can accurately measure your Vitamin D levels and assess whether you're getting enough sunlight. Our test uses liquid chromatography–tandem mass spectrometry (LC-MS/MS) – an advanced analytical method that ensures exceptional precision and sensitivity. Unlike standard Vitamin D tests, our method selectively measures only the active forms of Vitamin D and excludes inactive epimeric forms, which have a slightly different structure and do not share the same biological activity. Failing to exclude these epimers can lead to overestimated Vitamin D levels, especially in groups such as infants and pregnant women, where the proportion of epimers is naturally higher. If you require testing specifically for epimeric Vitamin D forms, please contact us for more information.

Specifications of *Bloodo™ Vitamin D Test*

Parameter	Specification	
	Vitamin D3	Vitamin D2
Method principle	Liquid chromatography with tandem mass spectrometry for Vitamin D3 and D2 detection.	
Suitability	For human dried blood spot analysis.	
Linearity	≥ 0.99	≥ 0.99
Error margin	29 %	30 %
Limit of detection, whole human blood	8.8 nmol/l	8.3 nmol/l
Limit of quantification, whole human blood	14.6 nmol/l	13.9 nmol/l
Linear range, whole human blood	12.5 – 249.6 nmol/l	12.1 – 242.3 nmol/l
Shipping conditions (DBS test kits)	Ship at ambient temperature.	
Storage conditions (DBS test kits)	Store at ambient temperature.	
Shipping conditions (DBS Test Cards, blood sample placed)	Ship at ambient temperature Post immediately. At the latest, on the following day. Sample stability limit – 29 days. Maximum sample shipping duration – 22 days.	
Home sampling conditions	Take the blood sample at a convenient time for you. Do not take supplements for 24 h before test. Dry the DBS Test Card fully for 3 hours.	

Key testing data of the *Bloodo™ Vitamin D Test*

1. Calibration and linearity

Vitamin D3 and D2 is measured according to a calibration curve received against the added Vitamin D3 and D2 standards in real human blood (Figure 1 and Figure 2). A typical calibration curve is conducted using Vitamin D3 and D2 concentration range of 12.5 – 249.6 nmol/l and 12.1 – 242.3 nmol/l (accounting for dilution), enabling the determination of Vitamin D3 and D2 in DBS samples within this range. Acceptance criteria for error margin is < 60 %. Error margin of the *Bloodo™ Vitamin D Test* for Vitamin D3 is 29 % and for Vitamin D2 – 30 %.

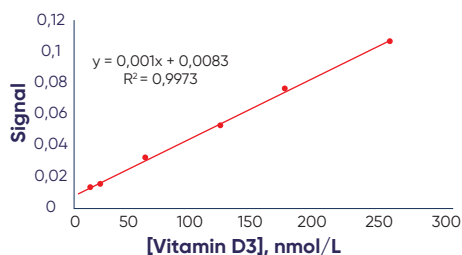


Figure 1. Typical calibration curve of *Bloodo™* Vitamin D Test. Vitamin D3 concentrations are given accounting for a dilution and shows total 25-hydroxyvitamin D3 metabolite concentration in human blood.

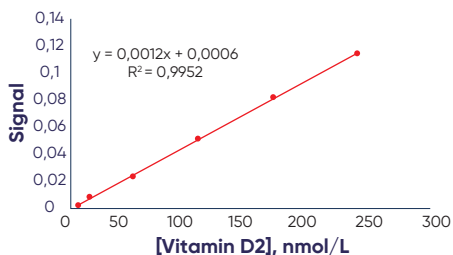


Figure 2. Typical calibration curve of *Bloodo™* Vitamin D Test. Vitamin D2 concentrations are given accounting for a dilution and shows total 25-hydroxyvitamin D2 metabolite concentration in human blood.

2. *Bloodo™* Vitamin D test limit of detection and limit of quantification

Limit of detection (LOD) is the lowest amount of analyte which can be detected, but not necessarily quantitated. The limit of quantification (LOQ) is the lowest concentration of an analyte in a sample that can be quantitatively determined with suitable precision and accuracy. LOD and LOQ for the vitamin D3 and D2 were determined from the linearity experiment results. The limit of detection and limit of quantification of the *Bloodo™* Vitamin D Test for vitamin D3 is 8.2 nmol/l and 14.6 nmol/l, respectively, for DBS samples (whole human blood). The limit of detection and limit of quantification of the *Bloodo™* Vitamin D Test for vitamin D2 is 8.3 nmol/l and 13.9 nmol/l, respectively, for DBS samples (whole human blood).

3. Specificity

Specificity is the ability of the assay to unequivocally assess the analyte in the presence of components that may be expected to be present. The characteristics of an LC-MS method impart remarkably high specificity for Vitamin D determination. High-performance liquid chromatography offers fast and effective initial separation of Vitamin D metabolites from the potentially interfering matrix components, reducing the complexity of analyte presentation to the mass spectrometer. Subsequently, mass spectrometry provides discriminatory power by filtering, detecting and quantifying molecules based on their parent ion distinctive mass-to-charge ratios (m/z). Moreover, tandem mass spectrometry (MS/MS), employing techniques like multiple reaction monitoring (MRM), enables fragmentation of the parent ion of Vitamin D metabolites and monitoring of highly specific product ions at different collision energies, further enhancing the method's specificity and minimizing the risk of misidentification. This methodology also selectively separates inactive C-3 epimeric forms of 25-hydroxyvitamin D, improving analytical accuracy and reducing the risk of overestimation, particularly relevant for infants, in whom epimer concentrations are higher.

NAD LEVEL TEST

Control your energy and longevity



Healthy
aging



Energy



Cognitive
function

- Track cellular energy
- Monitor aging process
- Evaluate mental clarity
- Check physical performance
- Analysis of NADH and NAD⁺

Order your
test here:



What is NAD?

Nicotinamide adenine dinucleotide (NAD) is a critical coenzyme in the human body. It exists in two interconvertible forms: the oxidized form NAD^+ and its reduced counterpart NADH. NAD is indispensable for energy production, powering muscle activity, respiration, and more. Beyond energy production, research shows that NAD also plays key roles in gene expression and DNA repair. As we age, cellular metabolism slows and DNA damage accumulates, increasing the demand for NAD and contributing to its gradual decline. Low NAD levels have been linked with higher risks of disease and cognitive impairment, and may manifest as symptoms such as fatigue and low energy.

Why test NAD?

NAD levels can decrease due to factors such as health conditions, lifestyle, and other individual differences, so testing can be beneficial for everyone. We strongly recommend testing if a person is using or planning to use NAD-boosting interventions. Currently, there are many different ways to raise NAD levels, including NR and NMN supplements, NAD injections, and IV drips. However, there are limited studies available on their actual long-term benefits and potential side effects. Without testing your baseline, supplementing with NAD is like shooting an arrow in the dark. Some supplements may not increase NAD levels at all, while others could raise them too quickly or excessively. The NAD test from dried blood is a safe and convenient way to evaluate your NAD status.

What kind of NAD do we measure?

NAD is present both inside and outside human cells. The balance between NAD^+ and NADH across different cellular compartments remains an active area of research, but studies suggest that overall NAD status can provide valuable insight into metabolic and cellular health. Our test measures total NAD in both forms and across all components of the blood. In other words, we measure the complete profile.

How does our test work?

A single droplet of blood is all it takes – the rest is up to us. Measuring NAD is not easy, since human cells like NAD and consume it constantly, so it is difficult to make an accurate analysis. Thousands of hours of scientific research were placed into the development of our NAD test. In order to detect the NAD, we have considered every variable, every chemical and every process – from blood sample collection on a special filter paper card to complex chemical reactions. The final result is a highly accurate, reproducible and selective test validated between two independent laboratories.

Specifications of *Bloodo™* NAD Test

Parameter	Specification
Method principle	Enzymatic cycling assay for total NAD (sum of NAD ⁺ and NADH) detection.
Suitability	For human dried blood spot analysis.
Relative standard deviation (RSD)	6.6 %
Limit of detection (LOD)	0.23 µmol/L
Limit of quantification (LOQ)	0.40 µmol/L
Linear range, whole human blood	5–60 µmol/L
Shipping conditions (DBS Test Kits)	Ship at ambient temperature.
Storage conditions (DBS Test Kits)	Store at ambient temperature.
Shipping conditions (DBS Test Cards, blood sample placed)	Ship at ambient temperature. Post immediately. At the latest, on the following day. Sample stability limit – 21 days. Maximum sample shipping duration – 14 days.
Home sampling conditions	Take blood sample as a first thing in the morning before having any food or drinks. Dry the DBS Test Card fully for 3 hours.

Key testing data of the *Bloodo™* NAD Test

1. Calibration and Error margin

NAD is measured according to a calibration curve received against the NAD standard (Figure 1). A typical calibration curve is conducted using NAD concentration range of 0–60 µmol/L (accounting for dilution), enabling the determination of NAD in DBS samples within this range. The LOD of the *Bloodo™* NAD Test is 0.23 µmol/L and LOQ is 0.40 µmol/L for DBS samples (whole human blood). Multiple DBS samples were analyzed and the relative standard deviation was calculated. The data yielded a 6.6% relative standard deviation.

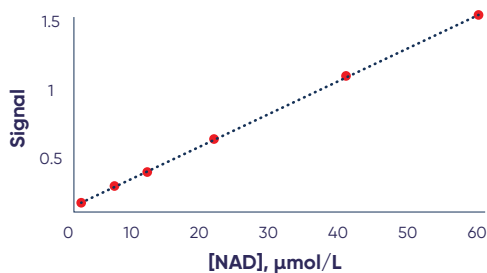


Figure 1. Typical calibration curve of NAD standard.

2. Bloodo™ NAD Test validation data

The Bloodo™ NAD test was validated by comparing the assay results between liquid and DBS samples ($R^2=0.94$) (Figure 2). The comparison demonstrated excellent correlation, indicating that the DBS samples align closely with the liquid blood samples. To validate that only NAD is being assayed, it was selectively removed from samples using NAD-specific enzymes (Figure 3). Normal samples showed typical NAD levels, while samples with NAD removed had significantly lower levels, confirming the Bloodo™ NAD Test's specificity.

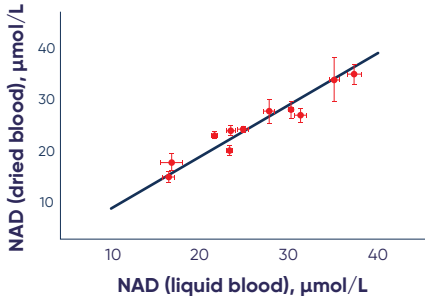


Figure 2. Comparison of testing results between liquid and dried blood samples using Bloodo™ NAD Test

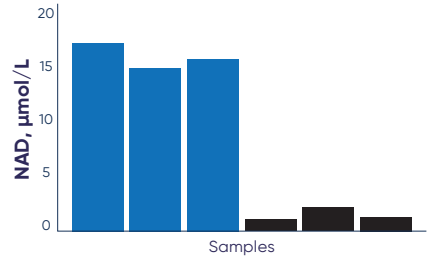


Figure 3. Testing NAD from typical DBS samples (blue) and samples after biochemical removal of NAD (black).

3. NAD supplement effectiveness monitoring using Bloodo™ NAD Test

Bloodo™ NAD Test was used to measure NAD levels from DBS samples taken from volunteers who were using NAD supplementation. Initially, all volunteers provided DBS samples for 5 days without any supplementation in order to establish baseline NAD levels in their blood. All samples were collected first thing in the morning before any food or drink intake. From the 5th day, three volunteers took 1 g per day of nicotinamide mononucleotide (NMN) in the morning, while the fourth volunteer, serving as a control, took 1 g per day of nicotinamide. The results are given in Figure 4 and show that Bloodo™ NAD Test measures NAD level changes in the blood after using NAD supplements.

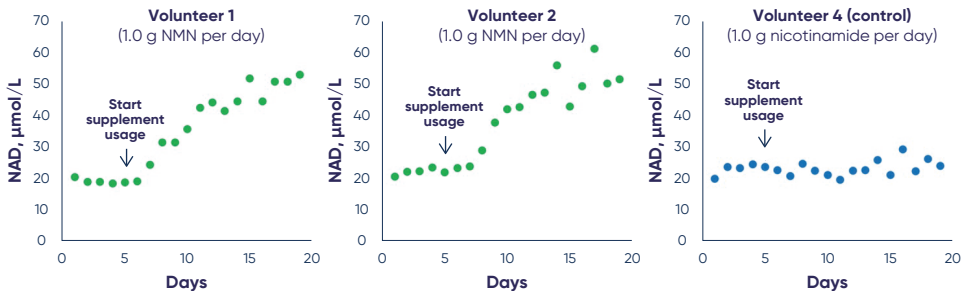


Figure 4. Bloodo™ NAD Test results of testing using DBS samples from volunteers taking NAD supplements.

HbA1c TEST

Control your long-term glucose levels



Diabetes
control



Diabetes
risk



Blood sugar
control


Order your
test here:



- Detect prediabetes and diabetes early
- Track your long-term blood sugar levels
- Optimize diabetes management
- Reduce risk of complications

What is HbA1c?

Hemoglobin A1c (HbA1c), or glycated hemoglobin, measures the percentage of glucose bound to hemoglobin in red blood cells. Hemoglobin is a protein found in red blood cells that carries oxygen throughout your body. This is an important blood test that gives a good indication whether a person is developing Type 2 diabetes or how well the diabetes is controlled.

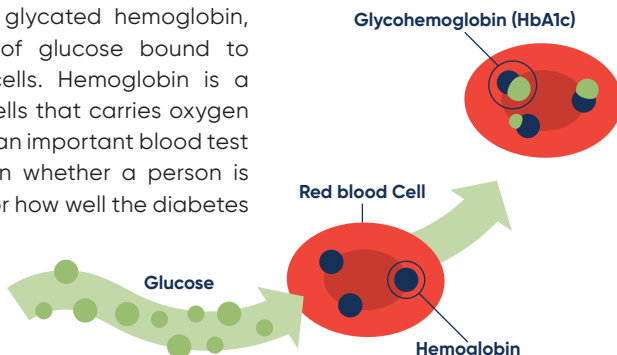

Why test HbA1c?

A high HbA1c means you have too much sugar in your blood. This means you're more likely to develop diabetes complications, such as serious problems with your eyes, feet or other conditions. Red blood cells have a life span of around 2-3 months, which is why we should monitor it quarterly. You can reduce your risk of health complications by knowing your HbA1c level and what you can do to lower it.

HbA1c testing applications:

- To identify individuals at increased risk for developing diabetes (pre-diabetic).
- Detection and treatment of pre-diabetes may reduce or eliminate the risk of developing Type 2 diabetes and related complications.
- Diagnosis of Type 2 diabetes.
- Improve treatments by maintaining glycemic control.

How does our test work?

The HbA1c test determines the concentration of glycated hemoglobin in the blood. The test is performed in a licensed clinical laboratory using the ion-exchange high-performance liquid chromatography method (HPLC) which is considered the gold standard for HbA1c analysis. HPLC method has the highest precision and accuracy among other methods used for HbA1c quantification. Two drops of blood is all it takes – the rest is up to us:

- No need to visit a hospital or a doctor.
- Easy and fast to perform.
- Minimally invasive, a simple finger prick is all it takes.
- Accurate and reliable results will be sent directly to you – take control of your own health!

Specifications of *Bloodo™ HbA1c Test*

Parameter	Specification
Method name	Ion-exchange high-performance liquid chromatography (HPLC) – gold standard for HbA1c testing. Test results are provided in NGSP % and IFCC mmol/mol units.
Suitability	For human dried blood spot analysis.
Repeatability	% CV 0.00-1.23 (NGSP %)
Accuracy	R ² =0,98
Linearity	From 3,9-18,8% (19-182 mmol/mol) with a maximum measured difference of ±0,1% (or ±1,0 mmol/mol) in this interval.
Shipping conditions (DBS Test Kits)	Ship at ambient temperature.
Storage conditions (DBS Test Kits)	Store at ambient temperature.
Shipping conditions (DBS Test Cards)	Ship at ambient temperature. Post immediately. At the latest, on the following day. Sample stability limit – 22 days. Maximum sample shipping duration – 14 days.
Home sampling conditions	Take the blood sample at a convenient time for you. Dry the DBS Test Card fully for 3 hours.

Key testing data of the *Bloodo™ HbA1c Test*

How are HbA1c results interpreted?

HbA1c values serve as pivotal indicators for the diagnosis of type 2 diabetes and the management of glucose control. The suspected diagnoses associated with HbA1c values, along with related recommendations, are detailed in Table 1.

HbA _{1c}		Suspected diagnosis	Recommendations
%	mmol/mol		
≥6,5	≥48	Diabetic	Immediately consult a doctor. Repeat the test every 3 months.
5,7-6,4	39-47	Pre-Diabetic	Increased risk of developing type 2 diabetes. Contact your doctor for more detailed examination, treatment and / or lifestyle correction. Repeat the test after 3 months.
<5,7	<39	Non-Diabetic	The risk of developing type 2 diabetes is low. In the absence of risk factors, repeat the test after 1 year.

Table 1. HbA1c result interpretation guide

Characteristics of the HbA1c test

1. Accuracy

HbA1c measured in capillary dried blood spot (CB DBS, x axis) samples were compared to three different sample matrixes to assess the accuracy of the technique.

To test each correlation, at the time of CB DBS collection, subsets of different donors provided other sample types:

- Whole blood EDTA (WB)
- Whole blood DBS (WB DBS)
- Capillary blood EDTA (CB)

Different sample matrix results are plotted against appropriate CB DBS sample results in Figure 1. R^2 values of each correlation are provided.

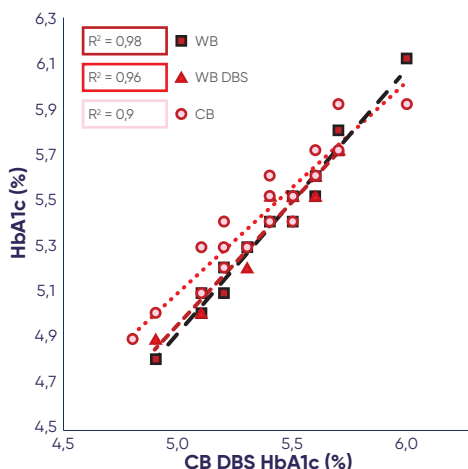


Figure 1. The accuracy of HbA1c results

2. Stability

For 22 days, HbA1c remains stable (%CV < 2 %) on the dried blood spot (DBS) card (Figure 2).

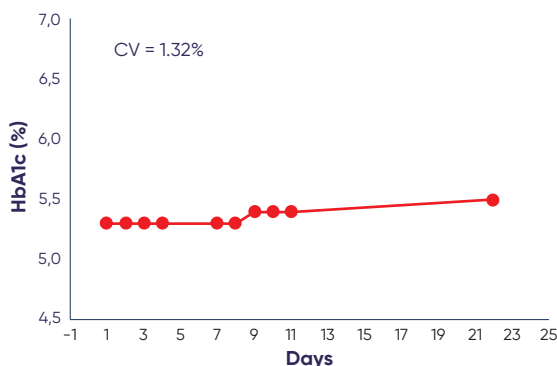


Figure 2. Stability of HbA1c on DBS cards

3. Repeatability

Between-run repeatability results were obtained for a series of samples by analyzing each sample 10 times in the same day. Coefficient of variation did not exceed 1.23% for all samples.

TESTOSTERONE+ TEST

Understand your inner fuel



Hormonal
balance



Energy &
resilience



Stress &
recovery



Libido &
well-being

- Assess testosterone status, hormonal balance, and bioavailability;
- Support energy, physical performance, recovery, and libido;
- Evaluate the impact of lifestyle, training, and hormone interventions;
- Detailed analysis of key hormones, including total testosterone, SHBG, calculated free testosterone, and complementary biomarkers.

What is Testosterone?

Testosterone is a key sex hormone present in both men and women, with naturally higher levels in males. It plays a crucial role in maintaining muscle mass, strength, energy levels, libido, and overall physical well-being. Beyond its physical effects, it also contributes to mood regulation, motivation, and cognitive function, supporting focus, confidence, and mental resilience.

In the bloodstream, testosterone exists in two main forms: bound (to sex hormone-binding globulin (SHBG) and albumin) and free (biologically active), which together determine its availability and biological effects in the body

Why test Testosterone?

Testosterone testing provides insight into androgen status, hormonal balance, and the bioavailability of testosterone in the body. Since only a small fraction of testosterone is biologically active, assessing total testosterone together with SHBG enables a more accurate estimation of free testosterone and its physiological effects.

Testosterone levels can be influenced by a range of physiological and pathological factors, including aging, metabolic health, endocrine function, and hormone-producing tissues such as the adrenal glands and gonads. Both low and elevated testosterone levels, as well as altered binding dynamics, may be associated with changes in energy, physical performance, body composition, libido, recovery, mood, and overall well-being.

How does our test work?

This test evaluates testosterone levels and their bioavailability using clinically validated, high-performance analytical methods. By integrating direct hormone measurement with advanced calculation models, it provides a clinically meaningful evaluation of androgen status and broader hormonal balance from a dried blood spot (DBS) sample. Total testosterone is quantified using liquid chromatography–tandem mass spectrometry (LC-MS/MS), the gold standard for steroid hormone analysis, ensuring exceptional specificity and sensitivity. SHBG is measured using a high-sensitivity enzyme-linked immunosorbent assay (ELISA) adapted for DBS samples. Free testosterone is then calculated using the Vermeulen equation, incorporating a standard albumin constant to estimate the biologically active hormone fraction.

Additional hormones, including DHEAS, cortisol, cortisone, androstenedione, and 17-OH progesterone, are measured using LC-MS/MS, enabling accurate assessment of adrenal function and overall hormonal balance.

What does the Testosterone+ Test measure?

Key biomarkers:

- **Total testosterone**
The overall amount of testosterone circulating in the bloodstream, including both bound and free forms.
- **SHBG (Sex Hormone-Binding Globulin)**
The primary transport protein that binds testosterone and regulates its availability in the body; higher SHBG levels are associated with lower free testosterone, while lower SHBG increases the free fraction.
- **Free testosterone (calculated)**
The biologically active fraction of testosterone, estimated from total testosterone and SHBG using the Vermeulen equation, a clinically validated model assuming standard albumin levels.

Additional biomarkers:


- **DHEAS**
Dehydroepiandrosterone sulfate is a hormone produced mainly by the adrenal glands and serves as a precursor to testosterone and estrogens, reflecting adrenal function, hormonal balance, and overall vitality.
- **Cortisol**
The body's primary stress hormone, regulating energy, metabolism, and the response to physical and psychological stress, while also playing an important role in modulating inflammation and maintaining physiological balance.
- **Cortisone**
The inactive form of cortisol, helping to regulate and protect the body from excessive cortisol activity, reflecting the balance between active and inactive stress hormones.
- **Androstenedione**
A naturally occurring androgen produced by the adrenal glands and gonads. It acts as a key precursor to testosterone and estrogens, indirectly contributing to the regulation of muscle mass, libido, reproductive health, and bone density.
- **17-OH Progesterone**
A steroid hormone and a key intermediate in the production of cortisol and androgens. It reflects the efficiency of hormone synthesis and supports the assessment of adrenal function.

Functional ratios

In addition to individual biomarkers, functional ratios provide deeper insight into hormonal balance and physiological regulation:

- **Cortisol / Cortisone**
Reflects the balance between active and inactive cortisol, providing insight into how cortisol is regulated and inactivated. This ratio helps assess the body's ability to control cortisol activity and maintain hormonal balance.
- **Cortisol / DHEAS**
Reflects the balance between stress response (cortisol) and anabolic or recovery-related processes (DHEAS). It provides insight into long-term adrenal function, resilience, and the body's ability to adapt to physical and psychological stress.
- **Testosterone / Androstenedione**
Reflects the efficiency of androgen conversion from precursor hormones to active testosterone. This ratio helps evaluate how effectively the body produces and utilizes testosterone within steroid hormone pathways.

Functional ratios reveal how hormones interact, not just their individual levels

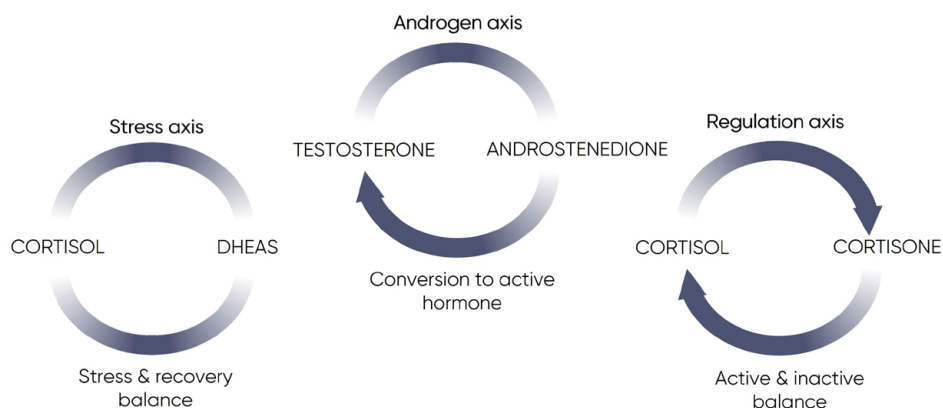

Specifications of *Bloodo™ Testosterone+ Test*

Parameter	Specification						
	Total Testosterone	SHBG	DHEAS	Cortisol	Cortisone	Androstenedione	17-OH Progesterone
Method principle	LC-MS/MS	ELISA	LC-MS/MS	LC-MS/MS	LC-MS/MS	LC-MS/MS	LC-MS/MS
Suitability	For human dried blood spot analysis.						
Linearity	≥ 0.99	≥ 0.99	≥ 0.99	≥ 0.99	≥ 0.99	≥ 0.99	≥ 0.99
Limit of detection, whole human blood	0.08 nmol/l	0.36 nmol/l	0.05 µmol/l	2.44 nmol/l	0.91 nmol/l	0.12 nmol/l	0.14 nmol/l
Limit of quantification, whole human blood	0.13 nmol/l	0.59 nmol/l	0.09 µmol/l	4.06 nmol/l	1.52 nmol/l	0.21 nmol/l	0.23 nmol/l
Linear range, whole human blood	0.16–41.6 nmol/l	0–295 nmol/l	0.15–190 µmol/l	8.64–1104 nmol/l	2.16–277 nmol/l	0.33–41.9 nmol/l	0.36–45.5 nmol/l
Precision (Inter-assay)	11.85 %	14.96 %	14.91 %	12.73 %	12.76 %	15.27 %	14.54 %
Low concentration level							
High concentration level	10.51 %	13.48 %	11.65 %	13.45 %	15.36 %	13.18 %	12.40 %
Shipping conditions (DBS test kits)	Ship at ambient temperature.						
Storage conditions (DBS test kits)	Store at ambient temperature.						
Shipping conditions (DBS Test Cards, blood sample placed)	Ship at ambient temperature Post immediately. At the latest, on the following day. Sample stability limit – 29 days. Maximum sample shipping duration – 22 days.						
Home sampling conditions	Take the blood sample in the morning (ideally between 7:00–10:00 AM, or within 1–2 hours after waking) to reflect peak testosterone levels. Avoid food for 8–12 hours before sampling. Avoid alcohol consumption, smoking, and intense physical activity for at least 24 hours prior to sampling. Avoid supplements that may affect hormone levels for 24 hours before testing. Dry the DBS Test Card fully for 3 hours.						

Technical details

1. Specificity

Specificity is ensured by LC-MS/MS through chromatographic separation and selective detection of each analyte based on its unique retention time and MRM transitions. Potential interferences from structurally related compounds and matrix components were evaluated and found not to affect analyte quantification. Stable isotope-labeled internal standards were used to correct for matrix effects and ensure accurate identification and quantification.

2. Trueness

Method trueness was evaluated across the working range and showed good agreement between measured and expected values. Recovery criteria (70–110%) were met for all analytes. At low and high concentration levels, recovery ranged from 80.9% to 103.2%, confirming accurate quantification across clinically relevant concentrations.

3. Precision

Precision was evaluated in terms of repeatability and within-laboratory reproducibility at low and high concentration levels. Repeatability ranged from 4.9% to 14.6% CV, while reproducibility ranged from 10.5% to 15.4% CV. All results were within the predefined acceptance criteria ($\leq 25\%$), confirming reliable and reproducible method performance.

4. Stability

Steroid hormones remained stable for up to 4 weeks (29 days), with measured concentrations ranging approximately from 85.1% to 104.8%, depending on the analyte and concentration level. The protein analyte SHBG also remained stable over the same period, with concentrations ranging from 80.2% to 109.3% across different concentration levels, all within predefined acceptance criteria.

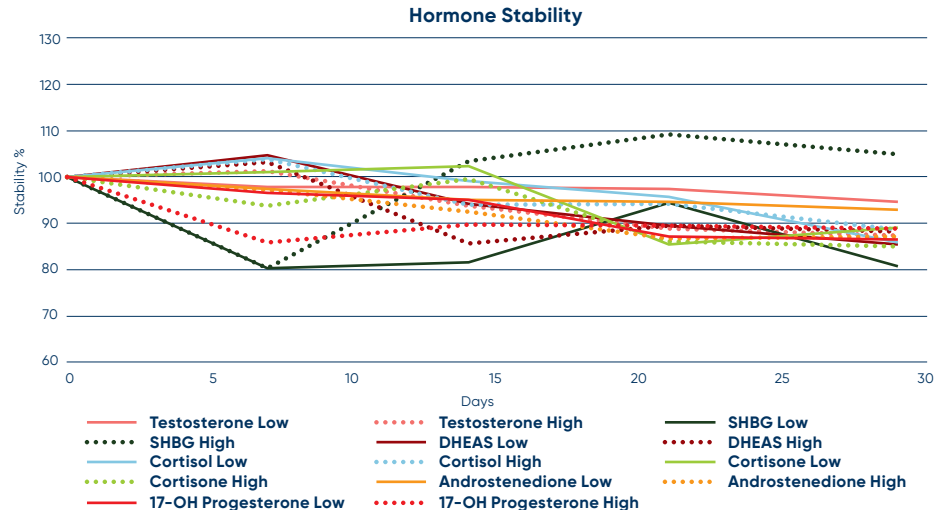


Figure 1. Hormone stability over 29 days at room temperature at low and high concentrations.



OMEGA RATIO TEST

Uncover your full essential fatty acids profile



Prick

>



Drip

>



Send



Design and components were carefully chosen to minimize environmental impact.

SAMPLE ID:



DBUQWNSQ2449



NAD LEVEL TEST

Control your energy and longevity



Prick

>



Drip

>



Send



Design and components were carefully chosen to minimize environmental impact.

SAMPLE ID:



DBUQW3PJ2449

 bloodo

HbA1c TEST

Control your long-term
glucose levels



Design and components were carefully chosen
to minimise environmental impact.

SAMPLE ID:



DBUQWNPQ248



4 779054 720481

 bloodo

VITAMIN D TEST

Check your sunshine levels



Design and components were carefully chosen
to minimise environmental impact.



4

779054 720481



4

779054 720672

